



ACCESSORIES & CONSUMABLES

ACCESSORIES

LAST UPDATED: June 2009

Glass Vials 5ml

This 5ml vial is manufactured from purest clear glass and is fully autoclavable. (186 vials per tray)

Code	Pack Size
VIALGLASS C-32	186



Silicone Bung 20mm

Made from the finest pure silicone, these bungs are fully autoclavable and designed for use with our glass vials and both types of aluminium crimp top seals.

Code	Pack Size
BUNG C-35	1500



Siloprene Bung 20mm

IOL grade. Made from the finest pure siloprene, these bungs are fully autoclavable and designed for use with our glass vials and both types of aluminium crimp top seals.

Code	Pack Size
BUNG SILOPRENE	1000

Horseshoe Crimp Top Seals (Various Colours)

Industry standard crimp manufactured from finest aluminium and available in a variety of colours.

Code	Pack Size
CRMPRED C-33A	1000
CRMP SILVER C-33A	
CRMPGOLD C-33A	
CRMPWHITE C-33A	
CRMPBLUE C-33A	



Flip Tear Crimp Top Seals (Various Colours)

Likely to become the new tamper evident standard for CE marked products. Made from pure aluminium with a colour coded plastic top.

Code	Pack Size
CRMPFLIP WHITE	1000
CRMPFLIP BLUE	





ACCESSORIES & CONSUMABLES

Plastic Tamper Evident Vial GP - 152

A well designed plastic Gas Permeable lens storage case. It features a soft plastic claw and screw top lid with a tamper evident lock.

Code	Pack Size
------	-----------

VIALGP – 152 T / E	250
--------------------	-----



Double Opticase Blue / White GP 151

This small compact Rigid Gas Permeable lens storage case features a plastic claw and two coloured lids with embossed L and R logo.

Code	Pack Size
------	-----------

OPTI GP151 B/W	250
----------------	-----



Double Soft Lens Mailers

A high quality double mailer available in clear plastic, affording excellent protection for soft lenses. This mailer can also be used for larger Gas Permeable lenses up to approximately 16mm in diameter

Code	Pack Size
------	-----------

MAILDOUBLE SL	100 and 1000
---------------	--------------



GP Double Mailers

These high quality double mailers are available in clear and are designed for use with Gas Permeable lenses. Can be cut for use as single units. *NB: Green is no longer available*

Code	Pack Size
------	-----------

MAILDOUBLE GP	250 and 1000
---------------	--------------




Rubber Sucker

This high quality two part green sucker is designed for handling Gas Permeable lenses.

Code	Pack Size
------	-----------

SUCKGREEN CM1059	10
------------------	----

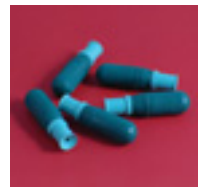

Scleral Suction Holder

This pure rubber suction holder is designed specifically for use with Scleral and other large diameter lenses.

Code	Pack Size
------	-----------

SUCKSCLERAL	3
-------------	---





ACCESSORIES & CONSUMABLES

Ultima Suction Holder

This high quality pure rubber suction holder is designed for Gas Permeable lenses. Supplied in individual storage cases.

Code	Pack Size
------	-----------

SUCKULTIMA	10
------------	----




White Plastic Tweezers

Large lens tweezers with super soft silicone tips for increased protection. This item can be over printed with your logo for a small extra charge, subject to minimum quantities. Also available at extra cost is a clear protection box.

Code	Pack Size
------	-----------

TWEEWHITE CH - 2074	250
---------------------	-----

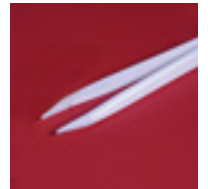


Yellow Flexible Tweezers 75mm

Small soft plastic yellow tweezers with a unique design that ensures constant tip pressure and therefore prevents lens damage. Ideal for retrieving soft lenses from glass vials.

Code	Pack Size
------	-----------

TWEEFLEX 75MM	100
---------------	-----

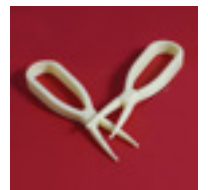

Shrink Wrap Sleeves

These clear shrink wrap sleeves are designed for use with our glass vials and feature a red tear away strip for easy tamper evident removal.

Code	Pack Size
------	-----------

SLEERED 37.5X45MM	1000
-------------------	------



Lens Catch Mat

This natural rubber mat is designed to be placed over any sink drainage hole to prevent the loss of your lenses.

Code	Pack Size
------	-----------

LENSMATS	1
----------	---





ACCESSORIES & CONSUMABLES

Steel Hand Crimper

Constructed from the finest materials, this high quality steel crimper is designed to give a lifetime's service. Available in versions suitable for Horseshoe or Flip tear crimps, it is the obvious choice for sealing vials.



Code	Pack Size
SHC 65 (Horseshoe)	1
FLIPCRIMPER (Flip Tear)	

Hydration Baskets

Available in various colours these baskets are used to hydrate soft contact and intraocular lenses. The pack sizes vary depending on the colour.



Code	Pack Size
HYDRcolour	Please Call Us

CONSUMABLES

Black High Melt Wax

This High Melt Wax is used for blocking Hydrophilic materials and is available in both 125g tubs and in packs of 10 sticks. Melting point 66°C.

Code	Pack Size
WAX BLACK HM 125G (Tubs)	125g
WAX BLACK HM STICK	10

Black Low Melt Wax

This Low Melt Wax is used for the blocking of Gas Permeable materials and is available in both 125g tubs and in packs of 10 sticks. Melting point 44 - 51°C.


Code	Pack Size
WAX BLACK LM 125G (Tubs)	125g
WAX BLACK LM STICK	10



C – 15 White Blocking Wax Sticks and Tubs

A Beeswax based blocking wax with a melting point between 62 and 64°C. This wax is used for blocking Hydrophilic materials and is available in packs of 20 sticks or in tubs of 500g.

Code	Pack Size
WAX WHITE C- 15 STICK	20
WAX WHITE C- 15 500G (Tubs)	500g

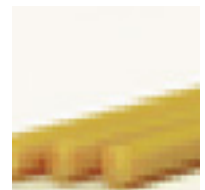



ACCESSORIES & CONSUMABLES

Kemdent Sticky Wax

This High Melt Wax is used for blocking Hydrophilic materials and is available in stick form in boxes of 500g. Melting point 55°C.



Code	Pack Size
WAX WHITE K 500G	500g

Contapol 1 - Soft Lens Polish

This is a high quality polishing compound suitable for all Hydrophilic materials, particularly those based on HEMA. It contains a 0.3 Micron uniformly dispersed alumina abrasive, suspended in an oil base, to produce excellent surfaces.

Code	Pack Size
CPOL 1	250g tubs

Contapol 2 - Soft Lens H W C Polish

This polish is specially formulated for use with high water content materials. It contains a 0.2 Micron abrasive suspended in a fine oil base, imparting excellent surface finishes. Available in 250g tubs.

Code	Pack Size
CPOL 2	250g tubs

Contapol 3 - Hard and GP Polish

This aqueous based polish contains a finely graded 0.2 Micron abrasive. Contapol 3 has been carefully formulated to eliminate solvent attack. A principle feature of the compound is the enhanced wettability imparted to materials polished with it. Suitable for all hard and Gas Permeable materials



Code	Pack Size
CPOL 3	250g tubs

Siflo – Hard and GP Polish

Siflo is a high quality alumina based polishing compound suitable for all hard and Gas Permeable contact lenses. The balanced viscosity ensures greater retention by the polishing cloth, reducing splatter and wastage. Siflo will not dry out during use, thus reducing the polishing time and effectively cooling the lens surface. Siflo rinses swiftly away in water leaving no residue.

Code	Pack Size
SIFL	250g tubs

PPOW .2 Micron

Polishing Powder for customers who wish to make up their own polishing compounds

Code	Pack Size
PPOW .2	500g and 10kg



ACCESSORIES & CONSUMABLES

Contasolve

This solution has been formulated for de-blocking and dissolving blocking wax on Gas Permeable and hard lenses.

Code	Pack Size
CSLV	1 litre cans



C – 18B White Brushed Cotton Cloth

This universal polishing cloth, widely used throughout the industry, combines excellent results with low cost.

Code	Pack Size
CLTHBRUSHED 18B	1 metre squares



C – 18C Beige Microfibre Cloth

The unique construction of the micro thin fibres is particularly suited to Gas Permeable lens polishing.

Code	Pack Size
CLTHBEIGE 18C	61cm squares 24X24 inches (610 x 610mm)



C – 18D Polishing Silk

This silk polishing cloth is designed for use on Gas Permeable and hard materials.

Code	Pack Size
CLTHSILK	1 metre squares



C – 19 White Delrin Polish Cups

A high quality plastic polishing cup designed for use with our Delrin cup sponges.

Code	Pack Size
DCUPWHITE C-19	10



C – 20 White Delrin Cup Sponges

A high quality foam sponge pads designed for use with our Delrin polishing cup.

Code	Pack Size
SPONWHITE C-20	100





ACCESSORIES & CONSUMABLES

C – 21 Roller Sponges

Designed for use with a variety of edge polishing machines.

Code	Pack Size
------	-----------

SPONYELLOW C-21	10
-----------------	----



CM - 1060 Foam Pads For Drum Tool

A high quality foam sponge pad designed for use in a Delrin drum polishing tool, for the modification or finishing of GP contact lenses.

Code	Pack Size
------	-----------

FPADWHITE 60	20
--------------	----



C – 22 Moulded Base Curve Polishing Cups

This high quality plastic moulded base curve polishing cup is available in black.

Code	Pack Size
------	-----------

MCUPBLACK 22	10
--------------	----

